




**SURE & EARLY**  
The most promising cycle days

Luxury Digital Ovation Test tells you 2 most fertile days



**What's special with Surcary™ Digital Ovation Test?**  
The Digital Ovation Test (D.O.T.) is accurate, taking away the confusion of interpreting lines.

• High chance of detecting your LH surge with 99% accuracy without compromising in terms of precision because of the LED display showing timing for the best and most accurate result.

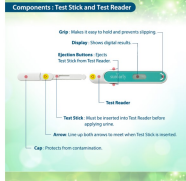

**How to Use**

**Step 1. Assemble the Test**

1. Take off the cap. (The test stick may be used immediately after removing the cap.)

2. Insert the stick. (Align the arrow on the Test Stick with the arrow on the Test Reader before inserting.)

• Assemble the test by inserting the Test Stick into the Test Reader.  
• Make sure the arrow on the Test Stick is facing the Test Reader's arrow.  
• Do not use the test if the arrow on the Test Stick is facing the Test Reader's arrow.

**Step 2. Test in your urine**

10. Remove the Absorbent Tip (which is used for Step 3).

3. Apply the urine stream directly on the Absorbent Tip for 3 seconds.

• No urinal test is needed for the Test Reader to work.

**Step 3. Wait**

• Wait for the Test Reader to show the test results.

• Do not touch the Test Reader or the Test Stick during the test.

• Do not use the Test Reader if the Test Stick is not inserted correctly.

**Step 4. Get your Results in about 3 Minutes**

**YES** (LH Surge) / **NO** (No LH Surge)

• YES: Indicates your LH surge has occurred. This means you are in your most fertile window. You should have sex now to increase your chances of becoming pregnant. You should have sex now to increase your chances of becoming pregnant.

• NO: Indicates you are not in your most fertile window. You should have sex later in your cycle to increase your chances of becoming pregnant.

• 99.9% accuracy in detecting your LH surge.

• Test message will remain on the Display for 7 minutes.

• Test message will remain on the Display for 7 minutes.

**Device management after use**

• **Disassembling the Test Stick:** Remove the cap from the Test Stick. Remove the Test Stick from the Test Reader. Dispose of the Test Stick in a safe and secure manner.

• **Disposing of the Test Reader:** Remove the Test Reader from the Test Reader. Dispose of the Test Reader in a safe and secure manner.

• **Disposing of the Test Reader:** Remove the Test Reader from the Test Reader. Dispose of the Test Reader in a safe and secure manner.

• **Disposing of the Test Reader:** Remove the Test Reader from the Test Reader. Dispose of the Test Reader in a safe and secure manner.

**Q&A**

**What does the 'No LH Surge' mean?**  
The 'No LH Surge' result indicates that you are not in your most fertile window. You should have sex later in your cycle to increase your chances of becoming pregnant.

**How accurate is the test?**  
The Digital Ovation Test has a 99.9% accuracy in detecting your LH surge.

**How long does the test take?**  
The test takes about 3 minutes to get your results.

**Can I use the test if I'm not in my most fertile window?**  
Yes, you can use the test at any time in your cycle to determine if you are in your most fertile window.

**Can I use the test if I'm not in my most fertile window?**  
Yes, you can use the test at any time in your cycle to determine if you are in your most fertile window.

**Can I use the test if I'm not in my most fertile window?**  
Yes, you can use the test at any time in your cycle to determine if you are in your most fertile window.